

Reagents

- PBS (Thermo Fisher Scientific, Cat#10010-23)
- Tris-HCl (Sigma, Cat#S7653)
- MgCl₂ (Sigma, Cat#63069)
- MgAc (Sigma, Cat#63052)
- NaCl (Sigma, Cat#S7653)
- KCl (Sigma, Cat#P9333)
- Protease Inhibitor (Roche, Cat#05056489001)
- RNase OUT (Invitrogen, Cat#10777-019)
- SUPERase IN (Invitrogen, Cat#AM2694)
- IGEPAL CA-630 (Sigma, Cat#I8896)
- Sucrose (Sigma, Cat#S7903)
- DTT (Sigma, Cat#D9779)
- Triton-X100 (Sigma, Cat#T9823)
- HEPES (Invitrogen, Cat#15630106)
- Spermidine (Sigma, Cat#85558)
- Digitonin (Sigma, Cat#D141)
- EDTA (Invitrogen, Cat#AM9261)
- dNTP (NEB, Cat#N0447S)
- Maxima Reverse H minus Reverse Transcriptase (Invitrogen, Cat#EP0751)
- T4 DNA Ligase (NEB, Cat#M0202L)
- NEBuffer 3.1 (NEB, Cat#B7203S)
- SDS (Invitrogen, Cat#15553-035)
- Protease K (NEB, Cat#P8107S)
- SPRI beads (Beckman coulter, Cat#B23319)
- Terminal Transferase (NEB, Cat#M0315S)
- KAPA HiFi HS (KAPA, Cat#KK2502)
- 10X Cutsmart buffer (NEB, Cat#M7204S)
- SbfI-HF (NEB, Cat#R3642)
- FokI (NEB, Cat#R0109S)
- NotI-HF (NEB, Cat#R3189)
- NEB Q5 DNA Polymerase (NEB, Cat#M0491)
- QIAquick PCR purification kit (QIAGEN, Cat#28104)
- KAPA qPCR library quantification kit (Roche)
- Barcodes plates and DNA primers, oligos (IDT)

The reagents in this section can be stored at -20 °C for up to 6 months

Barcode Plates (R02, R03 plates):

1. Distribute the R02 and R03 barcoded oligos to two 96-well plates, respectively as: 6 μL barcoded oligos (R02 or R03, 100 μM), 5.5 μL Linker (R02 or R03, 100 μM), 38.5 μL H_2O .
2. Seal the plates, heat for 5 mins at 95 $^{\circ}\text{C}$ and slowly cool down to 12 $^{\circ}\text{C}$ at the speed of 0.1 $^{\circ}\text{C}/\text{s}$.
3. Aliquot the each of the annealed barcodes plate into 4 working plates (each plate with 10 μL of annealed oligos in each well), sealed and stored at -20 $^{\circ}\text{C}$.

P5 adaptor mix:

1. Mix the oligo DNA as (a) P5 complex: 25 μL P5-FokI (100 μM), 25 μL P5c-NNDC-FokI (100 μM); and (b) P5H complex: 25 μL P5H-FokI (100 μM), 25 μL P5Hc-NNDC-FokI (100 μM).
2. Heat both stock mixture tubes for 5 mins at 95 $^{\circ}\text{C}$ and slowly cool down to 12 $^{\circ}\text{C}$ at the speed of 0.1 $^{\circ}\text{C}/\text{s}$.
3. Mixed P5 complex with P5H complex (on the ice) at the ratio of 1:3. Stock solution of P5-P5H becomes 50 μM , stored at -20 $^{\circ}\text{C}$.

Barcoded RT primers:

Mix the oligo (RNA barcode R01) in 12 corresponding 1.5 mL low-binding tubes: 12.5 μL RNA_RE (#01 to #12 in the 12 tubes, 100 μM), 12.5 μL RNA_NRE (#01 to #12 matched with RNA_RE, 100 μM), and 75 μL H_2O , stored at -20 $^{\circ}\text{C}$.

Barcoded Protein A-Tn5:

1. Mix the oligo (DNA barcode R01) in 12 PCR tubes: 25 μL DNA_RE (#01 to #12 in the 12 tubes, 100 μM), 25 μL pMENTs (100 μM , same in the 12 tubes).
2. Heat the mixtures for 5 mins at 95 $^{\circ}\text{C}$ and slowly cool down to 12 $^{\circ}\text{C}$ at the speed of 0.1 $^{\circ}\text{C}/\text{s}$.
3. Mix 1 μL of annealed transposome DNA with 6 μL of unloaded proteinA-Tn5 (0.5 mg/mL), briefly vortex and quickly spun down.
4. The mixtures were incubated at room temperature for 30 min then at 4 $^{\circ}\text{C}$ for an additional 10 min, stored at -20 $^{\circ}\text{C}$.

Tn5-AdaptorA:

1. Mix the oligo as: 25 μL Adaptor A (100 μM), 25 μL pMENTs (100 μM).
2. Heat the mixture for 5 mins at 95 $^{\circ}\text{C}$ and slowly cool down to 12 $^{\circ}\text{C}$ at the speed of 0.1 $^{\circ}\text{C}/\text{s}$.
3. Mix 1 μL of annealed transposome DNA with 6 μL of unloaded Tn5 (0.5 mg/mL), briefly vortex and quickly spun down.
4. The mixtures were incubated at room temperature for 30 min then at 4 $^{\circ}\text{C}$ for an additional 10 min. Dilute 10 X with dilution buffer (10 mM Tris-HCl pH7.5, 100 mM NaCl, 50% Glycerol, 1 mM DTT), stored at -20 $^{\circ}\text{C}$.

1 mM dCTP: 2 μ L dCTP (100 mM) + 198 μ L ultrapure H₂O.

The reagents in this section can be stored at 4 °C for up to 1 months

0.2% NIB (Nuclei isolation buffer):

1. Prepare the 10% IGEPAL CA-630 by 1:10 dilution with PBS.
2. Prepare the 0.2% NIB as: 200 μ L 10% IGEPAL CA-630, 0.5 g BSA, and add PBS to 10 mL.

2X TB (Tagmentation buffer):

1. Prepare the 2X TB as: 660 μ L Tris-Ac (pH 7.8, 1 M), 440 μ L K-Ac (3 M), 220 μ L Mg-Ac (1 M), 3.2 mL DMF, H₂O (5.48 mL).
2. Aliquot the 2X TB to 1.5 mL tubes and stored at 4 °C.
3. To avoid contamination, each experiment should use a new tube of 2X TB.

1M Spermidine: 145 mg in 1 mL ultrapure H₂O.

10% IGEPAL-CA630: 100 μ L + 900 μ L ultrapure H₂O.

2% Digitonin: 20 mg in 1 mL DMSO.

250 mM MgCl₂: 250 μ L MgCl₂ (1 M) + 750 μ L ultrapure H₂O.

The reagents in this section can be stored at RT for up to 1 months.

40.4 mM EDTA: 80.8 μ L EDTA (500 mM) + 919.2 μ L ultrapure H₂O.

Equipment

- Benchtop centrifuge (Eppendorf)
- ThermoMixer (Eppendorf)
- Swing bucket centrifuge (ThermoFisher)
- Qubit Fluorometric Quantification (ThermoFisher)
- Magnetic stands (96-well plate stand, PCR-tube stand)
- Fragment analyzer/Tapestation (Agilent)
- Micropipetts (0.5-20 μ L, 10-200 μ L, 100-1000 μ L)
- RNase-free tubes (300 μ L, 1.5 mL)
- Maximum recovery tubes (1.5 mL, Axygen, MCT150LC)
- 30 μ m Cell-Tric (Sysmex, Cat#04-004-2326)
- 96-well plates (Eppendorf, 0030603303)

Procedure

A. Day 1: nuclei preparation and antibody staining**A-1 Freshly prepare the following buffer:**

1. Prepare the Douncing Buffer (if processing mouse brain tissues) per sample as: 0.5 mL sucrose (1 M), 25 μ L KCl (2 M), 10 μ L $MgCl_2$ (1 M), 20 μ L Tris-HCl (1 M), 2 μ L DTT (1 M), 40 μ L Protease Inhibitor (50 X), 50 μ L SUPERase IN, 25 μ L RNase OUT and 1,328 μ L H_2O
2. Prepare the Complete Buffer as: 40 μ L HEPES (pH 7.5, 1 M), 75 μ L NaCl (4 M), 1 μ L spermidine (1 M), 40 μ L Protease Inhibitor (50 X), 50 μ L SUPERase IN, 25 μ L (RNase OUT), 2 μ L IGEPAL-CA630 (10%), 10 μ L digitonin (2%), 8 μ L EDTA (500 mM) and 1,749 μ L H_2O .
3. Prepare the 0.2% NIB-RPI (RNase and Protease Inhibitor) as: 1 mL 0.2% NIB, 20 μ L Protease Inhibitor (50X), 25 μ L SUPERase IN and 12.5 μ L RNaseOUT. Additional 0.2% NIB-RPI may be needed if there are multiple samples.

A-2 Nuclei preparation:I. Cultured cells:

1. Disassociated the cells according to culture protocol.
2. Spin-down cells at 300 g for 5 mins, resuspended the cells in 1 mL (5 million cells) 0.2% NIB-RPI. Sit on ice for 5 min.
3. Count the nuclei using cell counter.

II. Mouse brain tissues:

1. Prechill any tubes or tools, set the centrifuge to 4 °C.
2. Add 10 μ L 10% Triton-X100 in to the douncer (1 mL), add 1 mL Douncing Buffer.
3. Transfer the tissue dissection into douncer.
4. Use loose pestle 5-10 times gently followed by tight pestle for 15-20 times.
5. Filter with 30 μ m Cell-Tric. Transfer to 1.5 mL Maximum Recovery tubes. Spin down at 1000 g for 10 min.
6. Wash with 1 mL Douncing Buffer, spin down again at 1,000 g for 10min.
7. Resuspended the cells in 1 mL (5 million cells) 0.2% NIB-RPI. Sit on ice for 5 min.
8. Count the nuclei using cell counter.

A-3 Immuno-staining:

1. Aliquot 300 k nuclei into 12 Maximum Recovery tubes (total around 3.6 million nuclei); Species-mixing are recommended in at least one tube to estimate the potential collision rate.
2. Spin down the nuclei at 1,000 g for 10 min at 4°C and resuspended in 50 μ L Complete Buffer.
3. Add 2 μ g antibodies to each tube and incubate the mixture with rotation at 4 °C overnight.

B. Day2: multi-omics tagging and nuclei barcoding**B-1 Freshly prepare the following buffer:**

1. Prepare the Complete Buffer in 2 mL or appropriate tube as: 40 μ L HEPES (pH 7.5, 1 M), 75 μ L NaCl (4 M), 1 μ L spermidine (1 M), 40 μ L Protease Inhibitor (50 X), 50 μ L SUPERase IN, 25 μ L (RNase OUT), 2 μ L IGEPAL-CA630 (10%), 10 μ L digitonin (2%), 8 μ L EDTA (500 mM) and 1,749 μ L H₂O.
2. Prepare the Med Buffer #1 in 1.5 mL tube as: 20 μ L HEPES (pH 7.5, 1 M), 75 μ L NaCl (4 M), 0.5 μ L spermidine (1 M), 20 μ L Protease Inhibitor (50 X), 25 μ L SUPERase IN, 12.5 μ L (RNase OUT), 1 μ L IGEPAL-CA630 (10%), 5 μ L digitonin (2%), 4 μ L EDTA (500 mM) and 845 μ L H₂O.
3. Prepare the Med Buffer #2 in 5 mL or other appropriate tube as: 60 μ L HEPES (pH 7.5, 1 M), 225 μ L NaCl (4 M), 1.5 μ L spermidine (1 M), 60 μ L Protease Inhibitor (50 X), 75 μ L SUPERase IN, 37.5 μ L (RNase OUT), 3 μ L IGEPAL-CA630 (10%), 15 μ L digitonin (2%) and 2,532 μ L H₂O.

B-2 Nuclei washing and tagmentation:

1. Spin down nuclei at 600 g for 10 mins at 4°C, resuspended with 50 μ L Complete Buffer.
2. Repeat step 1 for 1 additional time.
3. Spin-down at 600 g for 10 mins at 4°C again, resuspended in 50 μ L Med Buffer #1.
4. Add 1 μ L pA-Tn5 (#01-#12) for each tube, rotate at room temperature for 60 min.
5. Spin-down nuclei at 300 g for 10 mins at 4°C, resuspended in 50 μ L Med Buffer #2.
6. Repeat step 5 for 2 additional times.
7. Add 2 μ L 250 mM MgCl₂, incubate in at ThermoMixer, 550 r.p.m., 37 °C, for 1 hr. During this time, pre-melt the barcoded RT primers on ice and mixed well.
8. Add 16.5 μ L 40.4 mM EDTA. Spin-down at 1,000 g for 10 mins at 4 °C. Discard the supernatant. Proceed to B-3 RT reaction immediately.

B-3 Reverse transcription (RT) reaction:

1. Distribute 4 μ L each of RT primers mix (#01-#12) into 12 PCR tubes.
2. Prepare the RT Mixture as: 52.8 μ L 5X RT Buffer, 60 μ L H₂O, 52.8 μ L PBS, 13.2 μ L dNTP, 3.3 μ L SUPERase IN and 1.65 μ L RNase OUT.
3. Resuspended the nuclei in 14 μ L RT mix and transfer to each corresponding RT tube with the same 1st round barcode in tagmentation step (R01).
4. Add 2 μ L of Maxima H- Reverse Transcriptase into each tube. Perform the following program in a thermocycler: step1, 50 °C * 10 mins; step2, 8 °C * 12 s, 15 °C * 45 s, 20 °C * 45 s, 30 °C * 45 s, 42 °C * 2 mins, 50 °C * 5 mins, repeat step 2 for additional 2 times; step 3, 50 °C * 10 mins; step 4, hold at 12 °C.
5. Combine nuclei from 12 tubes into a Maximum Recovery tube (ON ICE!!!), add 4.8 μ L 5% Triton X-100 and gently mix by pipette up and down.

6. Spin-down the nuclei at 1,000 g for 10 mins at 4 °C, discard the supernatant, resuspended in 1 mL 1X NEBuffer 3.1. Proceed to B-4 ligation-based barcoding immediately.

B-4 Ligation-based barcoding:

1. Prepare the of ligation mixture in 15 mL tubes as: 2,260 µL H₂O, 500 µL 10X T4 DNA Ligase Buffer, 50 µL BSA (20 mg/mL), 100 µL 10X NEBuffer 3.1. Two tubes of this mixture are needed.
2. Prepare the R02 Blocking Solution as: 264 µL R02 Blocker (100 µM), 250 µL 10X T4 DNA Ligase Buffer and 486 µL H₂O.
3. Prepare the R03 Termination Solution as: 264 µL R03 Quencher (100 µM), 500 µL EDTA (500 mM) and 236 µL H₂O.
4. Add the 1 mL nuclei (in 1X NEBuffer 3.1) into the 1st tube of Ligation Mixture. Add 100 µL T4 DNA Ligase into the mixture and mixed well.
5. Distribute the nuclei mixture into barcode plate R02: 40 µL nuclei mixture for each well. Seal the plate and incubate in a ThermoMixer at 37 °C, 300 r.p.m., for 30 mins.
6. Use multichannel pipette to add 10 µL of R02 Blocking Solution to each well and gentle mixed. Seal the plate. Continue the incubation for additional 30 mins in a ThermoMixer at 37 °C, 300 r.p.m..
7. Combine all the nuclei in a 15 mL tube (pre-washed with cold 0.5% BSA). Spin-down at 1,000 g, 10 °C for 10 mins. Discard the supernatant, resuspended in 1 mL 1X NEBuffer 3.1.
8. Add the 1 mL nuclei into the 2nd tube of Ligation Mixture. Add 100 µL T4 DNA Ligase into the mixture and mixed well.
9. Distribute the nuclei mixture into barcode plate R03: 40 µL nuclei mixture for each well. Seal the plate and incubate in a ThermoMixer at 37 °C, 300 r.p.m., for 30 mins.
10. Use multichannel pipette to add 10 µL of R03 Termination Solution to each well and gentle mixed. No further incubation is needed at this step.
11. Combine all the nuclei in a 15 mL tube (pre-washed with 0.5% BSA). Spin-down at 1,000 g, 4 °C for 10 mins. Discard the supernatant.
12. Wash the nuclei once with cold PBS. Spin-down at 1,000g 10 °C for 10 mins and resuspended in 200 µL - 1 mL cold PBS (optimal concentration 1,000 cell/µL). Proceed to lysis and DNA Cleanup.

B-5 Lysis and DNA Cleanup:

1. Count the nuclei, divide 2k-4k nuclei (2.5k recommended) into sub-libraries (e.g. 96-well plate).

Pause Point: Aliquoted nuclei can be stored at -80 °C for up to 6 months.

2. Prepare the lysis mixture as below: 3 µL 10% SDS, 3 µL Protease K, 3 µL NaCl (4M) and add PBS to total of 30 µL (including the volume of barcoded nuclei in PBS) for one sub-library. Make appropriate amount for the number of sub-libraries needed.

3. Add the lysis mix into nuclei. Seal the plate, incubate in a ThermoMixer at 55 °C, 850 r.p.m. for at least 2hr or overnight.
4. Cool-down the mixture to RT. Cleanup the DNA use 1.0X (30 µL) SPRI beads and elute in 12.5 µL Ultrapure H₂O. Store at -20 °C. During washing steps of the cleanup, please remove residual SDS as much as possible.

Pause Point: Purified DNA can be stored at -20 °C or -80 °C for up to 6 months.

C. Day3: library construction

C-1 TdT Tailing and pre-amplification

1. Add 1 µL 10X TdT Buffer and 0.5 µL 1 mM dCTP into 12.5 µL purified DNA. Incubate at 95 °C for 5 mins. Quickly chill on ice for 5 mins.
2. Add 1 µL of TdT to the mixture. Incubate at 37 °C for 30 mins followed by inactivation at 65 °C for 10 mins.
3. Add the Anchor Mix (6 µL 5X KAPA Buffer, 0.6 µL 10 mM dNTP, 0.6 µL 10 µM Anchor-FokI-GH, 7.3 µL H₂O and 0.6 µL KAPA HiFi HS DNA polymerase) to the reaction mix directly, and mix by pipette up and down. Perform the program in a thermocycler: step1, 95 °C * 3 mins; step2, 95 °C * 15 s, 47 °C * 60s, 68 °C * 120 s, 47 °C * 60 s, 68 °C * 120 s; repeat step2 for additional 15 times; step3, 72 °C * 10 mins; step 4, hold at 12 °C.
4. Add the Preamp Mix (4 µL 5X KAPA Buffer, 0.5 µL 10 mM dNTP, 2 µL 10 µM PA-F, 2 µL 10 µM PA-R, 11 µL H₂O and 0.5 µL KAPA HiFi HS DNA polymerase) to the reaction mix directly and mix by pipette up and down. Perform the program in a thermocycler: step1, 98 °C * 3 mins; step2, 98 °C * 20 s, 62 °C * 20 s, 72 °C * 150 s; repeat step 2 for additional 9 times; step 3, 72 °C * 120 s; step4, hold at 12 °C.
5. Use 0.2X + 0.65 X (10 µL + 32.5 µL) of SPRI beads for double-size selection purification. Elute in 36 µL ultra-pure H₂O, use 1 µL for Qubit quantification. Recode the concentration x ng/µL. Depends on the cell type and numbers of cells per sub-library, in our hand typically x will in the range of 1 ng/µL to 30 ng/µL (see **Troubleshooting** for details).

Pause Point: Purified DNA can be stored at -20 °C or -80 °C for up to 6 months.

C-2 2nd adaptor tagging and digestion

1. Transfer 17 µL of the purified amplified product into two 200 µL PCR tubes for DNA and RNA, respectively.
2. For DNA part: add 2.5 µL 10X Cutsmart buffer, 1 µL SbfI-HF, 1 µL FokI, and 3.5 µL H₂O.
3. For RNA part: add 2 µL 10X Cutsmart buffer, 1 µL NotI-HF.
4. Incubate at 37 °C for 60 mins.
5. Use 1.25 X SPRI beads (31.3 µL for DNA, 25 µL for RNA). Elute both DNA and RNA in 10 µL H₂O.

Pause Point: Purified DNA can be stored at -20 °C or -80 °C for up to 6 months.

6. For DNA part: add 2 μ L 10X T4 DNA Ligase Buffer, 1.5 μ L P5 Adaptor Mix, 5 μ L H₂O and 1.5 μ L T4 DNA Ligase. Incubate with the ligation program (4 °C * 10 mins, 10 °C * 5 min, 16 °C * 15 min, 25 °C * 45 min).
7. For DNA part: Cleanup the ligation mixture with 1.25 X (25 μ L) SPRI beads and elute in 30 μ L H₂O.
8. For RNA part: add 10.5 μ L 2X TB, y μ L 1:10 diluted Tn5-AdaptorA ($y = x * 0.2$, *titration may be needed for different batch of Tn5*).
9. For RNA part: The tagmentation reaction was carried out 550 r.p.m., 37 °C for 30 min in a ThermoMixer followed by cleaned up using QIAquick PCR purification kit and eluted in 30 μ L 0.1X EB (QIAGEN).

Pause Point: Purified DNA can be stored at -20 °C or -80 °C for up to 6 months.

C-3 Indexing PCR and sequencing

1. For DNA libraries, prepare the PCR mix as: 30 μ L purified template, 2 μ L 10 μ M TruSeq i7 indices, 2 μ L 10 μ M Universal P5 primer, 1 μ L 10 mM dNTP, 10 μ L 5X NEB Q5 Buffer, 4 μ L H₂O and 1 μ L Q5 HiFi DNA polymerase. Run the program in a thermocycler: step1, 98 °C * 3 mins; step2, 98 °C * 10 s, 63 °C * 30 s, 72 °C * 60 s; repeat step2 for 8 cycles; step3, 72 °C * 60s; step4, hold at 12 °C.
2. For RNA libraries, prepare the PCR mix as: 30 μ L purified template, 2 μ L 10 μ M TruSeq i7 indices, 2 μ L 10 μ M NextEra N5 indices, 1 μ L 10 mM dNTP, 10 μ L 5X NEB Q5 Buffer, 4 μ L H₂O and 1 μ L Q5 HiFi DNA polymerase. Run the program in a thermocycler: step1, 72 °C * 5 mins, 98 °C * 30s; step2, 98 °C * 10 s, 63 °C * 30 s, 72 °C * 60 s; repeat step2 for 8 cycles; step3, 72 °C * 60s; step4, hold at 12 °C.
3. Use 0.5 μ L (dilute 1000X) to run a qPCR to estimate the additional cycles needed (typically 0-4) to reach 10 nM. Run the additional cycles in a thermocycler: step1, 98 °C * 3 mins; step2, 98 °C * 10 s, 63 °C * 30 s, 72 °C * 60 s; repeat step2 for desired cycles; step3, 72 °C * 60s; step4, hold at 12 °C.
4. Cleanup the libraries using 0.9 X (45 μ L) SPRI beads, elute in 15-30 μ L H₂O. Use qPCR to quantify libraries. Store at -20 °C.

Pause Point: Purified libraries can be stored at -20 °C or -80 °C for up to 6 months.

5. Use Agilent Fragment Analyzer to analysis the fragment distribution of libraries (**Fig. 1**). If there is a clear band at ~245 bp (**Fig. 1c**), repeat the size selection with 0.75X SPRI beads or other size selection method (see **Troubleshooting** for details).
6. Multiplex libraries and sequence with standard Illumina sequencing primers on commercial sequencing platforms. We have tested sequencing of Paired-Tag libraries on HiSeq 2500/4000, NextSeq 500 or NovaSeq 6000 platforms. Libraries should be loaded at recommended concentrations according to manufacturer's instructions. At least 50 and 100 sequencing cycles are needed for Read1 and Read2, respectively. For example: using PE 50 + 7 + 100 cycles (Read1 + Index 1 + Read2) on a NextSeq 500 platform with 150-cycle sequencing kits, or PE 100 + 7 + 100 cycles on a NovaSeq 6000 platform with 200-cycle sequencing kits.

D. Initial data processing and QC

1. All scripts for initial Paired-Tag data processing are available from GitHub (<https://github.com/cxzhu/Paired-Tag>). Download and perform the data pre-processing steps according to instructions on the GitHub repository.
2. Initial Paired-Tag data processing includes (a) extract barcode sequences from Read2, (b) assign barcode combinations to cellular barcodes references (assign barcode sequences to ID of 12 sample tubes and 2 rounds of 96 wells), (c) mapping the assigned reads to reference genome and (d) generating cell-to-features matrices for downstream analyses.
3. The following metrics during initial Paired-Tag data processing can be used for quality control. For step 2(a), typically >85% and >75% of DNA and RNA reads will have full ligated barcodes. For step 2(b), >85% of both DNA and RNA reads can uniquely assigned to one cellular barcode with no more than 1 mismatch. For step 2(c), typically >85% of assigned reads can be mapped to the reference genome; depending on which histone mark targeted, from 60% to >95% of assigned DNA reads can be mapped to the reference genome.
4. Cell-to-features matrices generated from initial data processing steps is a standard sparse matrix format which can be used for down-stream analyses with commonly used single-cell analysis software.

Troubleshooting

- After pre-amplification, typically ~30 – 1,000 ng DNA will be recovered (around 1 ng/ μ L to 30 ng/ μ L measured by Qubit, when pre-amplified products were eluted in 36 μ L). If the concentration is less than 0.2 ng/ μ L: (1) considering check the reagents used for pre-amplification steps (C1); (2) pre-amplification may be inhibited by residual SDS from B-5 step4, considering repeat this step (B-5 step4); (3) the quality of extracted nuclei may not be good enough for Paired-Tag experiment, considering perform a bulk ATAC-seq from extracted nuclei for QC. If the concentration is higher than 50 ng/ μ L, the amplified products are likely to be over-amplified barcode-adaptor dimers, please see the next paragraph for details.
- If only bands at ~245 bp are observed for DNA and RNA libraries, the libraries are likely failed as barcode-adaptor dimers without any insertions are the dominant component. Possible solutions include (1) considering optimization of tissue dissociation and nuclei preparation steps, (2) removing the supernatant as much as

possible in B-2 step9, B-3 step6, B-4 step7 and step11; (3) considering repeating the B-4 step12 once if enough nuclei (>100k) are recovered.

- If both bands of ~245 bp and 300-1000 bp exist for DNA or RNA libraries (**Fig. 1c**), indicating the presence of both libraries and barcode-adaptor dimers. (1) If a clear ~245 bp band exists, consider perform an additional 0.75X SPRI beads size selection (C-3 step4) (or other size selection method) to remove the primer-barcode dimers (~245 bp). (2) If there is no ~245 bp band but with significant fraction of fragments with length < 300 bp (e.g., **Fig. 1c**), consider reduce the amount of N5-Tn5 used in C-2 step8.

Time taken

Nuclei extraction (A-1 to A-2): 2 hr

Immuno-staining (A-3): overnight

Nuclei washing and tagmentation (B-2): 3.5 hr

Reverse transcription (B-3): 1.5 hr

Nuclei barcoding (B-4): 2 hr

Lysis and DNA cleanup (B-5): overnight or 3 hr

Preamplification (C-1): 4 hr

Adaptor tagging (C-2): 3 hr

Indexing PCR and QC (C-3): 3-5 hr

Anticipated results

The genomic DNA and cDNA libraries can be generated from desired number of nuclei, depending on how many sub-libraries were used for library construction. In our hands, the lengths of DNA fragments in genomic DNA libraries are typically between 300bp

- 700bp (average ~400bp), for those from cDNA libraries are mainly between 300 bp
- 1000 bp (average ~550bp) (**Fig. 1a, b**).

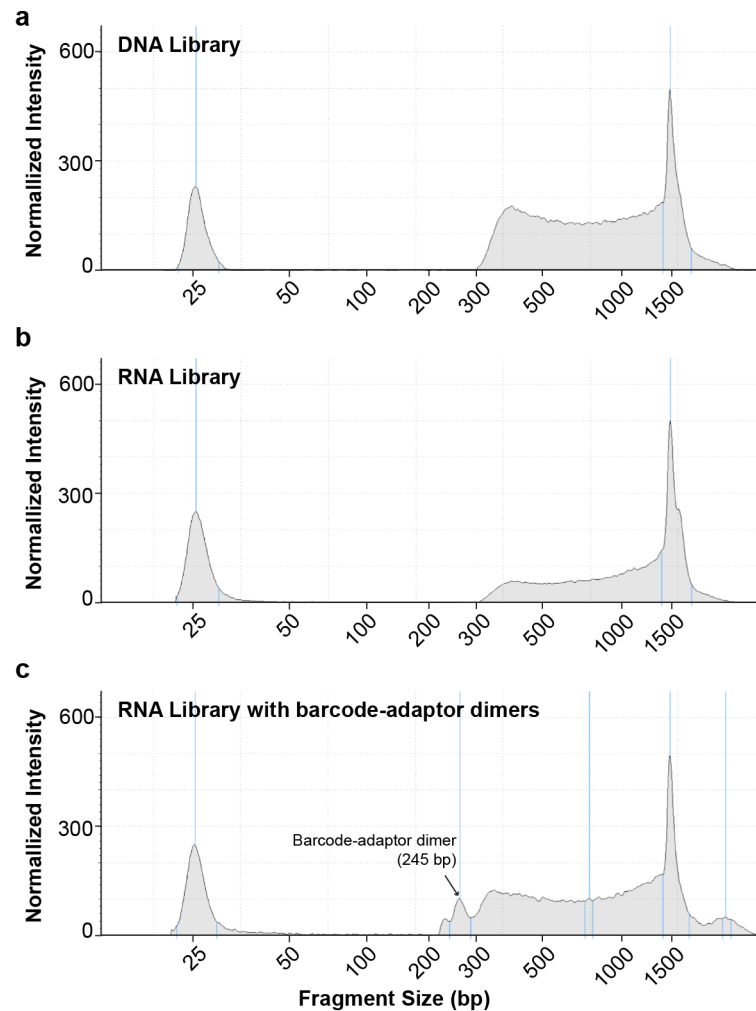


Fig.1 | Representative fragment size analysis results of Paired-Tag libraries.

Results of Agilent 2200 TapeStation system with High Sensitivity D1000 ScreenTape from (a) a Paired-Tag DNA library, (b) an RNA library and (c) a representative RNA library with barcode-adaptor dimers that need to be removed.